



HIGH INTENSITY
GAMMA SOURCE



UNIVERSITY OF VIRGINIA



TRIANGLE
UNIVERSITIES

HI γ S FROZEN SPIN TARGET

HiFrost Manual

Author:
Ryan DUVE

Supervisors:
Blaine NORUM
Don CRABB
Pil-Neyo SEO

Updated June 2, 2014

Preface

HiFrost is nuclear polarized target apparatus consisting of a dilution refrigerator, internal magnetic coil, microwave guide and NMR coil. External components of HiFrost include a polarizing magnet, microwave generating EIO, pump and vacuum system to run the dilution refrigerator, and the Q-Meter/Yale Card set up for running the NMR.

To polarize a target, chemically doped or irradiated beads are placed in the inner-most chamber, the fridge is assembled around the beads and the whole fridge is swung to align with the HIGS gamma beam. The beads are cooled with liquid helium and polarized with a large external magnet, with an EIO providing microwaves necessary for dynamic nuclear polarization.

To achieve frozen spin mode, the internal magnet is ramped on as the external magnet is ramped off, so eventually the external magnet can be powered down and removed while keeping the polarized beads in a steady magnetic field.



Figure 1: The HiFrost refrigerator during a test cooldown at UVa.

Contents

1	Theory of Operation	1
1.1	Nuclear Polarization	1
1.1.1	DNP	1
1.2	NMR	1
1.3	Frozen Spin	3
1.4	Dilution Refrigerator	3
2	List of Hardware	5
2.1	Fridge	5
2.1.1	Assembly Tools	5
2.1.2	Cryogenic Instruments	6
2.1.3	Magnet Leads	7
2.2	Target Material	7
2.3	Fridge Stand Mounting/Swinging Tools	7
2.4	Pump Station Area	8
2.4.1	^3He Pumps (and oil)	8
2.4.2	^4He Pumps	8
2.4.3	LN_2 Trap	8
2.4.4	Helium Gas Cylinders	8
2.5	Vault Pumps	9
2.5.1	OVC	9
2.5.2	IVC	9
2.6	Cryogenics	9
2.6.1	Safety Equipment	9
2.6.2	Transfer Lines	9
2.7	Vault Electronics Rack	10
2.7.1	Modules	10
2.7.2	Patch Panel	10
2.7.3	Thermistor System	10

2.8	NMR	11
2.9	Microwaves	12
3	Safety	13
3.1	Cryogenic Safety	13
3.1.1	Equipment	14
3.1.2	Cryogenic Explosions	14
3.1.3	Cryogenic Burns	15
3.2	Oxygen Depravation	17
3.3	High Voltage Safety	17
3.4	Gravity Safety	17
3.5	Magnetic Field Safety	19
3.6	Chemical Safety	19
3.6.1	Heavy elements	19
3.6.2	Isopropyl Alcohol	20
4	Practical Operation	21
4.1	Coldown	22
4.1.1	Prep	22
4.2	LHe Transfer Overview	22
4.2.1	Dewar Filling	22
4.2.2	Cool ^4He Section	29
4.2.3	Cool ^3He Section	30
4.2.4	Dilution	30
4.3	Polarization	30
4.4	NMR	30
4.5	Frozen Spin	30
5	Procedures	33
5.1	Evacuating Transfer lines	33
5.1.1	Cryofab TLs	33
5.1.2	LTL	35
5.2	Indium Extrusion	39
5.2.1	Background	39
5.2.2	Materials	40
5.2.3	Extrusion process	41
5.3	Making Indium Seals	42
5.4	Breaking Indium Seals	44
5.5	Weighing LHe Dewars	45
5.6	Magnet Lead Installation	45

5.7	Mounting Fridge	45
5.8	Dismounting Fridge	45
5.9	Cleaning Viton/Buna O-rings and Grooves	45
5.9.1	Fridge flanges	46
5.10	Changing Kenol Connectors	47
5.11	Installing Heating Tapes	47
5.12	Filling LN ₂ Trap	47
5.13	Check 250LD/500LD LHe Level	47
5.14	Attach Pressurization Line to LHe Dewars	47
5.15	Purging Helium Lines	48
5.16	Vacuum Rise Testing	48
5.17	Stripping Fridge	49
5.18	Assembling Fridge	49
5.18.1	MC	50
5.18.2	Microwave Guide Support	50
5.18.3	IVC	52
5.18.4	OVC	52
5.19	Removing/Replacing ³ He Baffles	53
6	Subsystems	55
6.1	Thermistor Switch and Display	55
Appendices		
Appendix A Checklist for cooldown		59
Appendix B SLPM Conversion		61
Appendix C Transfer Line Drawings		63

Chapter 1

Theory of Operation

1.1 Nuclear Polarization

In the presence of a magnetic field, spin-1/2 nuclei tend to align themselves along the axis of the field. The polarization of the ensemble of particles is defined by

$$P = \frac{\uparrow - \downarrow}{\uparrow + \downarrow},$$

where (\downarrow , \uparrow) is the population of nuclei with $m_z=(-1/2, 1/2)$.

1.1.1 DNP

Dynamic Nuclear Polarization (DNP) is a process of using microwave radiation to pump electron-proton pairs to higher energy (polarized) states. The relaxation time for the electron is much shorter than that of the proton, so the proton remains polarized while the electron is able to be paired with other protons for polarization.

The microwave frequency can either be the sum of the proton and electron splitting frequency or the difference, depending on whether the nucleons are to be aligned with or against the magnetic field.

1.2 NMR

Nuclear Magnetic Resonance (NMR) is the method we use to measure the polarization of the target material. Essentially, a coil in the immediate vicinity of the target acts as the inductor in an LCR circuit, the inductance



Figure 1.1: DNP diagram [CM97]. The white arrows are electron polarization and the black arrows are proton polarization. Microwaves at a frequency difference $\nu_e \pm \nu_p$ flip the spins of both particles, and the electron's superior lattice coupling ensures it will flip back before the proton.

of which changes as a function of the magnetic susceptibility of the nuclei in the coil. Since the magnetic susceptibility and polarization are related, the resonance frequency of the LCR circuit tell us information about the polarization of the target material[CG80].

1.3 Frozen Spin

The target is polarized in a 2.5 T field, called the polarizing field. The target will lose all polarization in an environment free of magnetic field. However, the detector cannot fit around the polarizing magnet, and the polarizing magnet would obstruct outgoing neutrons, anyway). To maintain polarization while the target is in the detector, an internal 0.5 T field is turned on inside the fridge and the polarizing magnetic is wound down and removed. The 0.5 field, called the holding field, sufficiently maintains polarization of nucleons while the detector is moved around the target and data is taken.

1.4 Dilution Refrigerator

Since the polarization goes like the inverse of temperature, colder environments make for better polarized targets. In our case, we use a dilution refrigerator, one that mixes the two isotopes ^3He and ^4He for cooling, to reach target temperatures below 0.1 K.

The dilution refrigerator principle relies on the splitting of a $^3\text{He}/^4\text{He}$ mixture into two distinct phases, a ^3He concentrated phase and a ^3He dilute phase. The area labeled “Two-phase region” in Figure 1.2 illustrates which mixtures, characterized by ^3He concentration, are inaccessible at which temperatures. Since two distinct phases in thermal contact are always in or striving for thermodynamic equilibrium, changing the concentration of the dilute phase (by pumping ^3He out of it) will cause atoms from the concentrated phase to cross the phase boundary to restore balance. Since the heat change of mixing (enthalpy difference between the dilute phase and concentrated phase) is positive, the ^3He crossing the phase boundary must absorb energy from the surrounding environment, which it does in the form of heat.[Hoc]

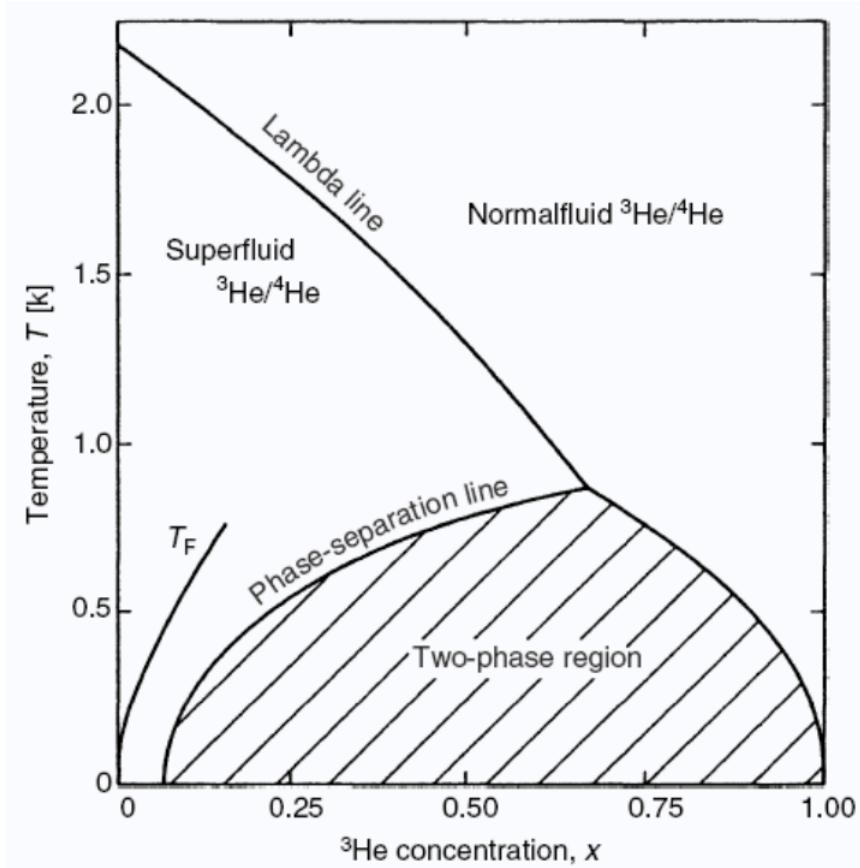


Figure 1.2: The famous diagram that shows the splitting of two distinct phases of ${}^3\text{He}-{}^4\text{He}$ mixture. [Bet89]

Chapter 2

List of Hardware

This chapter is an index of all the equipment used in the project. Part and/or serial numbers are included where applicable.

2.1 Fridge

2.1.1 Assembly Tools

OVC

- wire gripping pliers
- 5 mm allen key
- tee pilot pins
- 8x M6x25 triple flange machine screws

IVC

- solder station
- 2.5 mm allen key
- indium waste box
- wooden or hard plastic q-tips
- wire cutting pliers
- isopropyl alcohol
- 6x M3x12 machine screws, 316SS only

Microwave guide support

- flat head screwdriver
- fine tip squeeze-to-open tweezers

Mixing chamber

- mixing chamber customized open wrench
- 1.5 mm allen key
- indium waste box
- wooden or hard plastic q-tips
- wire cutting pliers
- isopropyl alcohol
- 6x M2x12 machine screws and nuts, 316SS only

2.1.2 Cryogenic Instruments

^4He sensors

- 2x *LakeshoreTM* DT-670 silicon diode
- *LakeshoreTM* PT-102 platinum RTD
- *LakeshoreTM* Quad-Lead cryogenic wire, 36 AWG
- *LakeshoreTM* Varnish, IMI-7031

^3He sensors

- *LakeshoreTM* DT-670 silicon diode
- *OmegaTM* KHLV-0502 Kapton heater
- *LakeshoreTM* RX-202A-AA ruthenium oxide sensor
- *LakeshoreTM* GR-200A-30 germanium diode
- 4-wire Matsushita resistor (model unknown)

2.1.3 Magnet Leads

- OFHC copper leads
- insulating mylar
- Kapton tape
- small flat head screwdriver
- minature flat head screws
- *Solid SealingTM* KT20678 copper conductor KF feedthrough

2.2 Target Material

- cryogenic storage LN₂ dewar
- open styrofoam tub, about 5-10 liters
- tweezers
- hemostats
- LN₂ transfer dewar
- cryogenic gloves

2.3 Fridge Stand Mounting/Swinging Tools

- 13 mm combination wrench
- long-shaft 3/16" ball driver
- 10 mm combination wrench
- 5 mm allen key
- fridge lever
- *80/20TM* 5/16"-18 nuts, bolts

2.4 Pump Station Area

- Alcatel or VIC leak checker
- KF-25 plastic blank off
- KF-25 flexible hose
- Adixen 40 HPX oil mist eliminator

2.4.1 ${}^3\text{He}$ Pumps (and oil)

- *PfeifferTM* WKP4000 Roots pump (P3)
- *AlcatelTM* MIV1000 Roots pump
- *AlcatelTM* MIV350 Roots pump (VP1)
- *AlcatelTM* 2063-H mechanical pump (A119)
- *AlcatelTM* ZM2004H mechanical pump (A114)

2.4.2 ${}^4\text{He}$ Pumps

- *AlcatelTM* RSV1002 Roots pump (A120)
- *AlcatelTM* RSV300 Roots pump (A100 or A113)
- *AlcatelTM* 2063-H mechanical pump (A119)

2.4.3 LN₂ Trap

- LN₂ transfer dewar
- cryogenic gloves
- heat gun

2.4.4 Helium Gas Cylinders

- 1-1/8" open wrench
- 1/2" fNPT to hose barb fitting for dewar pressurization

2.5 Vault Pumps

2.5.1 OVC

- *PfeifferTM* Hi Cube 80 Eco,PM S03 555 A turbo pump
- *KeyTM* flexible 316SS hose, KF25, 77"
- *PfeifferTM* PT R25 500 ActiveCold Cathode gauge
- OVC manifold assembly

2.5.2 IVC

- turbo pump
- *PfeifferTM* PCR 280, ActivePirani/Capacitance gauge
- *KeyTM* flexible 316SS hose, KF16, 78"
- IVC manifold assembly

2.6 Cryogenics

2.6.1 Safety Equipment

- cryogenic handling gloves
- closed toed shoes
- apron
- safety goggles
- ODH monitors

2.6.2 Transfer Lines

Cryofab TLs

- WTL
- UTL
- MTL
- *CryofabTM* OPE1/2-KF25 SV-9 Operator KF25 pumpout adapter

Cern LTL

- LTL
- LTL pump out adapter
- teflon tape
- flat head screwdriver
- LTL oring *McMasterTM* #9452k26 Buna-N AS568A Dash #114

2.7 Vault Electronics Rack

2.7.1 Modules

- 2x *LakeshoreTM* 218S readout
- AVS47 resistance bridge
- *CryomagneticsTM* Model 1200 S/H level readout
- *AgilentTM* E3646A power supply
- Holding coil power supply
- Polarizing magnet power supply

2.7.2 Patch Panel

- small flat head screw driver
- small philips head screw driver

2.7.3 Thermistor System

- *OmegaTM* #SA1-TH-44006-T EIO body thermistor
- *OmegaTM* #ON-909-44006 water flow thermistors
- *NewportTM* INFCH-011 thermistor meter/controller
- *MCCDAQTM* 1608G USB DAQ
- *OmegaTM* DP25TH thermistor readout

- *Gems Sensors*TM M103005 (product #02H2390 from *Newark*TM) flow meter display
- *Gems Sensors*TM 155421 water flow meter

2.8 NMR

1. computer with PDP installed
2. *NI*TM BNC-2090 rack mounted terminal block
3. water cooled housing
4. *Bay Voltex*TM HT-5000-AC-RC-VR water chiller
5. power supply
6. DIO panel
7. SMA wrenches
8. *Rhode and Schwartz*TM SMA100A function generator
9. cable from R&S to SMA
10. *NI*TM PCIe-GPIB 778930-01 computer card
11. *NI*TM 763061-02 cable
12. *NI*TM PCIe-6321 - 781044-01 computer card
13. *NI*TM 192061-02 cable
14. *NI*TM PCIe-6509 - 779976-01 computer card
15. *NI*TM 182762-02 cable
16. $\lambda/2$ cable
17. oscillator crystals
18. plastic screwdrivers
19. Q-Meters
20. Yale Cards
21. oscilloscope

2.9 Microwaves

- EIO
- Cober power supply
- *Phase MatrixTM* 578 microwave counter
- *HPTM* 432A power meter
- *PhasematrixTM* ACC093 mixer
- various microwave guides and attenuators
- 5/64" ball driver
- 3/32" ball driver
- water chiller
- *OmegaTM* SA1-TH-44006-40-T thermistor
- *OmegaTM* DP25-TH thermistor readout
- *Gems SensorsTM* 155421 flow sensor
- *Gems SensorsTM* M103005 M103 series flow readout

Chapter 3

Safety

This chapter is about personnel safety, not equipment safety. Learning to use HiFrost equipment properly without damaging anything is the aim this whole manual. Safety of lab users, however, is important enough to warrant an entire devoted chapter, even if many points are reiterated throughout the document.

Ways to die or become seriously injured working on HiFrost:

- Cryogenic explosions
- Cryogenic burns
- Oxygen deprivation
- Electric shock
- Falling
- Hot surfaces
- Chemical exposure

3.1 Cryogenic Safety

Both liquid nitrogen (LN_2 , 77 K) and liquid helium (LHe , 4 K) are used during HiFrost operation. The two primary cryogenic safety concerns

are over-pressurization (explosion) of cryogenic vessels and cryogenic burns. LN₂ and LHe are each at risk of both, and precautions are taken while handling either.

3.1.1 Equipment

Cryogenic gloves

Cryogenic gloves are lined internally with super-insulation and are long enough to cover past halfway between the wrist and elbow. HiFrost equipment includes large and medium sized sets of gloves, and the best fitting size should always be used. Typical activities requiring use of cryo gloves are filling the LN₂ trap, inserting and removing LHe transfer lines, filling LN₂ dewars for the target loading procedure, and making or manipulating target beads.

Goggles

Wearing goggles and/or a face shield while handling open LN₂ dewars is recommended. Goggles also protect users in the vicinity of the large helium gas plumes expected during LHe dewar filling procedures.

Closed toed shoes

The general FEL safety regulations preclude anyone from wearing open toed shoes in controlled areas. Still, there are areas outside the radiation perimeter where cryogens may be handled (e.g., the supply dewar just outside the bay doors), and closed toed shoes must be worn when handling LN₂ or LHe vessels.

3.1.2 Cryogenic Explosions

The principle hazard of handling cryogens is an explosion caused by an enclosed volume of liquid warming up without adequate pressure relief for the evaporated gas. For this reason, every container that holds cryogenic liquid has a pressure relief valve, and most large dewars have a non-configurable, one time use emergency burst disk.

An important exception to the pressure relief rule is the ³He circuit in the dilution refrigerator and pumping system. Due to the scarcity and cost of ³He, it is unacceptable to install pressure relief valves venting to atmosphere. Instead, a single internal pressure relief valve, which opens

at about 1.2 bar, is installed on the $^3\text{Hegas}$ rack and should open if the pressure in the circuit rises due to a plug. The pressure relief valve can only open if the valves leading to the fridge (via the condenser and still lines) are opened in the appropriate configuration. Failure to do this could lead to catastrophic damage to the pumping system or, worse yet, the refrigerator itself.

Three valves rule

The 500LD, 100LD and most supply helium dewars have three pathways for helium to escape in addition to the burst disk. The three pathways are:

1. the pressurization port, where external gas is applied to the dewar
2. a pressure relief valve, usually set around 2-4 PSI
3. the inlet/outlet pathways where transfer lines connect through the top of the dewar

In general, all three of these ports may have manual valves to prevent helium from flowing through them. For example, at the end of a 500LD fill, we close the outlet port on the filling dewar to halt the transfer, and the pressure relief valve is already closed to maintain liquid flow to the 500LD (see Section 4.2.1). The pressurization valve is always closed unless specifically venting the dewar or actively applying pressure to it. The three valves rule is

At least one of the three pathways on a cryogenic vessel must be open at all times.

If all three pathways are closed simultaneously, the radiative heat load (present in all dewars) will warm up the trapped liquid, increasing the pressure on the dewar until the burst disk breaks. If for whatever reason the emergency burst disk fails, the resulting pressure will lead to an enormous explosion, easily fatal to any personnel nearby [Mat06].

3.1.3 Cryogenic Burns

Cryogenic burns may happen by touching a liquid cryogen or cold gas plume, touching a bulk mass that was recently cooled by cryogenic liquid, or inhaling cold gas. Use cryogenic gloves whenever handling liquid cryogens or surfaces they have recently cooled (like transfer lines), and always wear closed toed shoes in the lab.

In the event of a cryogenic burn[Epl14]:

- If prehospital warming is attempted, options include placing the affected area in warm (not hot) water or warming it using body heat (eg, placing frostbitten fingers in the axillae).
- You can flushing the area with tepid water but, in order to avoid tissue damage, a forceful flow of water should NOT be used.
- Never use dry heat. Never apply direct heat.
- Do not rewarm frostbitten tissue if there is a possibility of re-freezing before reaching definitive care. This would result in worse tissue damage.
- Do not rub frostbitten areas in an attempt to rewarm them; this can cause further tissue damage.
- Remove any clothing or jewelry that may restrict circulation to injured area.
- Do not cover the area; leave injured area open to air.

To get help[Epl14], call Duke Employee Occupational Health and Wellness Hotline, available at any time, for assistance with a typical burn, 919-684-8115. A Duke Hospital Operator will answer and then will page the EOHW nurse to your number. During business hours you may also call EOHW main number, 919-681-3136, option 2 and ask for nurse. Call 911 if situation is life threatening. Provide the following information[Epl14]:

- Agent causing cryogenic burn
- Skin tone/color (white, pale, red, or blue).
- Presence of Blisters (yes or no)
- Location and area of injury
- Presence of open skin (yes or no)
- Current treatments

- Any additional past medical information, medications, allergies
- Confirm that the area of liquid nitrogen or hydrogen is secure and safe

3.2 Oxygen Depravation

Work with cryogenics automatically involves an oxygen depravation hazardous (ODH) environment, because the contents any liquid cryogen vessel are capable of displacing many times its volume of breathable air. Generally, liquid helium can displace 1 cubic meter of breathing air for each liter of liquid that is quickly boiled, meaning the 100 L HiFrost dewar can displace 100 cubic meters of air. Additionally, helium is colorless, odorless and tasteless, so it is often impossible for a worker to tell when they are not getting enough oxygen until the symptoms of oxygen deprivation begin to kick in.

The following information and the content of Table 3.1 is taken from Thomas Jefferson National Lab's ODH manual. [TJNL]

Health Effects of Reduced Oxygen Normal air is approximately 21% oxygen and 78% nitrogen. The remaining 1% is mostly argon. Health effects begin at an oxygen concentration of 17%. Oxygen monitors at Jefferson Lab are set to alarm at 19.5%. This advance warning should give ample time to escape the hazard area. The early health effects are difficult to detect so the oxygen monitors are relied upon to give early warning:

3.3 High Voltage Safety

3.4 Gravity Safety

No high places should be reached without the appropriate ladder during Hifrost operations or activities. Remember the following 5 rules of ladder safety[OO]:

1. Choose the right ladder for the job.
2. Inspect the ladder before you use it.
3. Set up the ladder with care.
4. Climb and descend ladders cautiously.
5. Use safe practices when working on a ladder.

Additional information ladder safety is provided by OSHA[SA].

Percent Oxygen	Health Effects
17	night vision reduced increased breathing volume accelerated heartbeat
16	dizziness reaction time for new tasks is doubled
15	poor judgement poor coordination abnormal fatigue upon exertion loss of muscle control
10-12	very fault judgement very poor muscular coordination loss of consciousness
8-10	nausea vomiting coma
<8	Permanent brain damage
<6	spasmodic breathing convulsive movements death in 5-8 minutes

Figure 3.1: JLab ODH manual's list of health effects of oxygen deprivation.

3.5 Magnetic Field Safety

Aside from being full of 77 K liquid nitrogen and 4 K liquid helium, the polarizing magnet poses a safety threat by virtue of its 2.5 T field. The 0.5 T frozen spin field is considerably weaker and, due to its localization inside the radiation shield of the fridge, less accessible to personnel and foreign objects.

The most likely danger is a ferromagnetic object (like a wrench) striking personnel as it is pulled towards the polarizing magnet[IC]. The best safety measure against this is making sure Hifrost workers and anyone else in the GV know when the magnet is on and do not keep loose ferromagnetic objects on their person.

Anyone with an implanted medical device, such as a pacemaker, should alert Hifrost officials before working on or near Hifrost equipment during a polarizing magnet cooldown. Under no circumstances should they get closer than the 10 gauss line while the magnet is energized without speaking with their physician and the director of DFELL[RWDL].

The EIO has a strong enough magnetic field to pull a ferromagnetic object in from a few inches away, so care should be taken to not allow any loose tools to smash into it.

3.6 Chemical Safety

3.6.1 Heavy elements

Indium

Indium is used for sealing two interfacing sets of flanges in the dilution refrigerator. Disassembling the refrigerator necessarily entails scraping indium off these flanges, and assembling the fridge requires cutting and setting the seal from a roll of indium wire.

While no cases of indium poisoning by oral consumption have been recorded, it remains a heavy element and toxic to humans if it enters the blood stream. As a precaution, wearing latex gloves is recommended while handling indium, and hands should be washed immediately after, especially before taking a lunch break or leaving for the day.

Lead

Lead is a toxic metal often used in the laboratory due to its density and availability. Inorganic lead is not readily absorbed through the skin, but

once it enters the body (through ingestion or inhalation) it is carried by the bloodstream to the “bone, teeth, liver, lungs, kidneys, brain and spleen” in high concentrations [Sta98]. Children and pregnant women are particularly at risk to damage from lead poisoning [EPA13]

Lead bricks are used at the University of Virginia PTGroup lab to counter the weight of the fridge on the stand. Lead bricks are also present at HIGS around the beam line and UTR, although they are not usually found in the Vault. Work gloves and hard-toed shoes are recommended when carrying or lifting lead bricks. Any skin that was exposed to lead bricks or lead dust should be washed immediately after performing duties requiring lead exposure.

3.6.2 Isopropyl Alcohol

Isopropyl alcohol is used to clean the indium and KF-oring surfaces on the fridge and in the pumping system. It is generally safe to be exposed to, but if there is a large quantity in an open container (like a bath for soaking vacuum parts) make sure the area is well ventilated and anyone working nearby knows about it. Symptoms of isopropyl alcohol inhalation are dizziness, drowsiness and headache, and may cause unconsciousness [Air].

Isopropyl alcohol has a flash point (the lowest temperature it emits ignitable fumes) of 53° F, and care should be taken not to expose alcohol bottles to open flames.

Chapter 4

Practical Operation

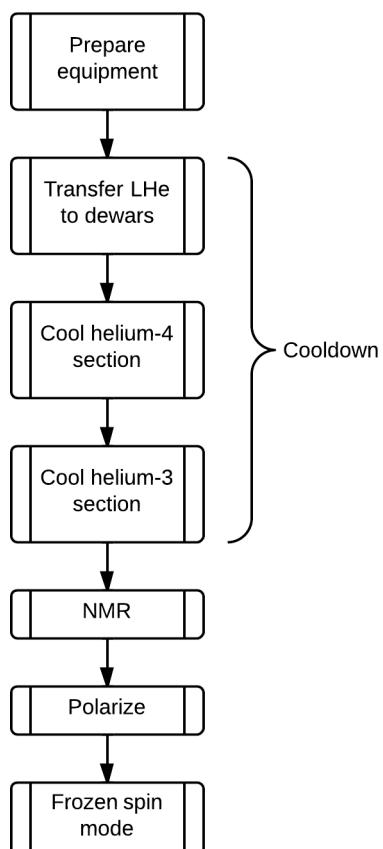


Figure 4.1: Practical steps to polarize target material with HiFrost.

This chapter assumes the reader is familiar with everything safety related. Read Chapter 3 before attempting anything here.

4.1 Cooldown

4.1.1 Prep

Go down the checklist found in Appendix A, reporting any issues to the senior cooldown scientist.

This step is not optional: the rest of this operational procedure assumes all prep work is successfully completed and equipment bugs are worked out. See walkthroughs for specifics tasks in Chapter 5.

4.2 LHe Transfer Overview

Figure 4.2 is a schematic of the LHe transfer system. More detailed figures of the dewars are found later in the chapter, and complete technical drawings of the TLs are found in Appendix C.

4.2.1 Dewar Filling

The 500LD and 100LD are initially filled before cooling the fridge. The magnet is filled either before or after dilution.

The 250LD has an LHe outlet port, a pressurization/exhaust port and the pressure relief port, shown in Figure 4.3. Each port has an accompanying valve.

The 500LD has an inlet port (for LHe filling and the level probe), an outlet port, a pressurization port, a vacuum pump out port and a pressure relief port, as shown in Figure 4.4. Only the pressurization, vacuum pump out and pressure relief ports have accompanying valves.

500LD Fill

1. Weigh a full 250LD according to the procedure in Section 5.5 before breaking the seal on it.
2. Check 250LD level using the procedure in Section 5.13.
3. Place correct Goddard fittings on the 250LD and/or the UTL.

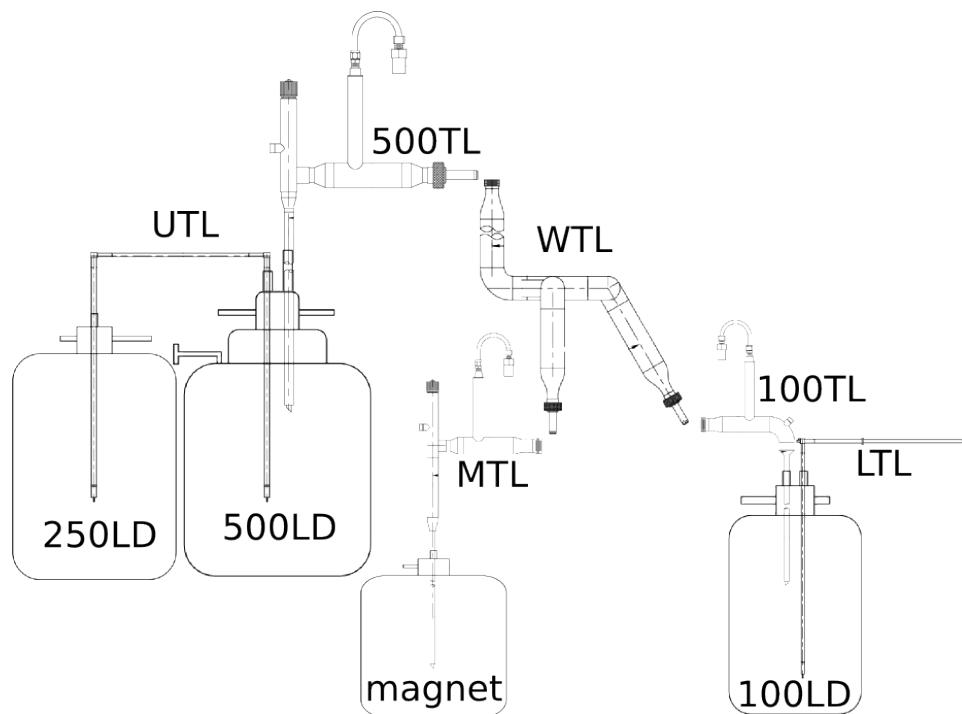


Figure 4.2: Schematic overview of the LHe transfer and storage system.

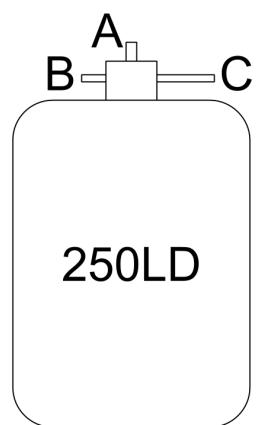


Figure 4.3: An example of 250LD ports: A) LHe outlet port, B) pressure relief valve, C) pressurization/exhaust port.

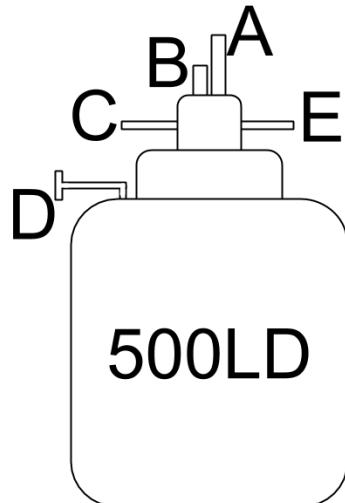


Figure 4.4: The 500LD ports.

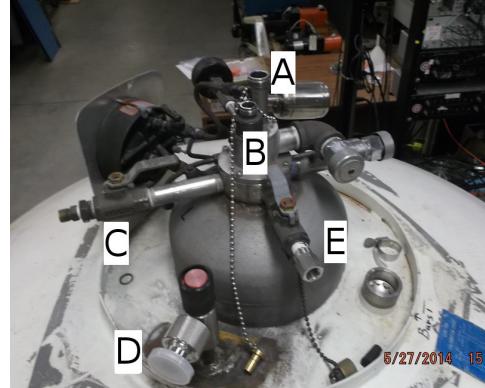


Figure 4.5: Photo of 500LD: A) LHe inlet port, B) LHe outlet port, C) pressurization port, D) vacuum pump out port, E) pressure relief valve.

4. Blow out UTL with helium gas by fitting the universal hose over it and putting about 2 PSI on the helium cylinder.
5. [initial 500LD fill only] Purge the 500LD with helium gas by hooking up the pressurization line, closing the main flow valve, opening the pressure release valve and putting 5 pounds of pressure on the 500LD for about 10 minutes (Figure 4.7). Remove the pressurization line.
6. Open main flow valve on 500LD to vent pressure (if any), then open pressurization valve making sure the pressurization port is not pointed at anyone.
7. Release the 250LD pressure from the main flow valve and close the pressure relief valve. Slowly lower the correct end of the UTL down into the 250LD about 30 cm or until helium gas is blowing out the 500LD end, as shown in the left side of Figure 4.8. Then, lower the 500LD end in as far as it can go, while continuing to slowly lower the 250LD side at a rate of about 2 cm/s. A cold plume should be coming out of the 500LD. Check that no gas is escaping through the Goddard fittings on either dewar (if it is, replace fitting o-rings).
8. Hook up the pressurization line to the 250LD pressurization port using the procedure in Section 5.14 and open the pressurization valve as

shown in the right side of Figure 4.8. There should now be 2 PSI on the 250LD.

9. Periodically measure the 250LD level using the procedure in Section 5.13. When it is empty, stop flow from the pressurization cylinder regulator and disconnect the pressurization line.
10. Simultaneously raise both sides of the UTL out of the dewars and hang it on the wall to warm up.
11. Bung the 500LD inlet port, close the 500LD pressurization valve, and leave open the pressure relief valve.
12. On the 250LD, open the pressure release valve and close the main flow and pressurization valves.
13. Recover the Goddard fittings from the 250LD and bring it back to the weigh area. Record the final weight on the log sheet.

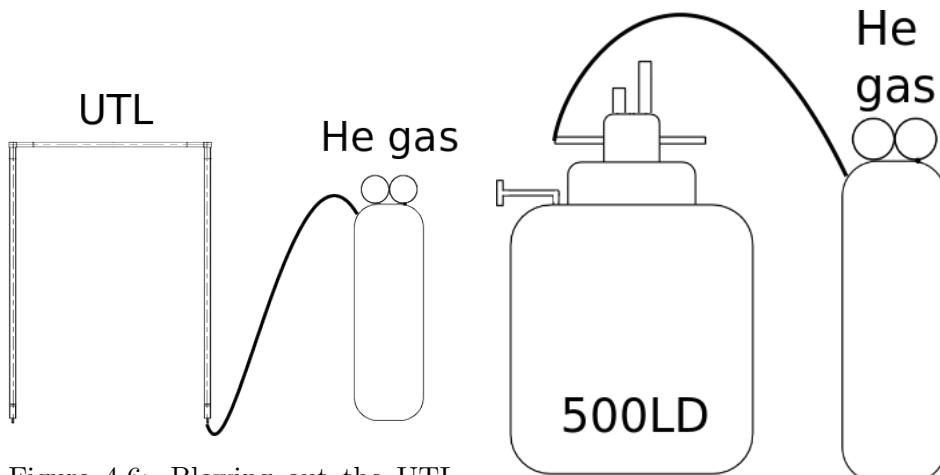


Figure 4.6: Blowing out the UTL with He gas (Step 4).

Figure 4.7: Initial 500LD fill, purging the 500LD with He gas (Step 5).

Initial 100LD Fill

The 100LD has an entrance port, a main flow port, an exhaust port and a pressure relief port. The exhaust and pressure relief ports have accompanying valves.

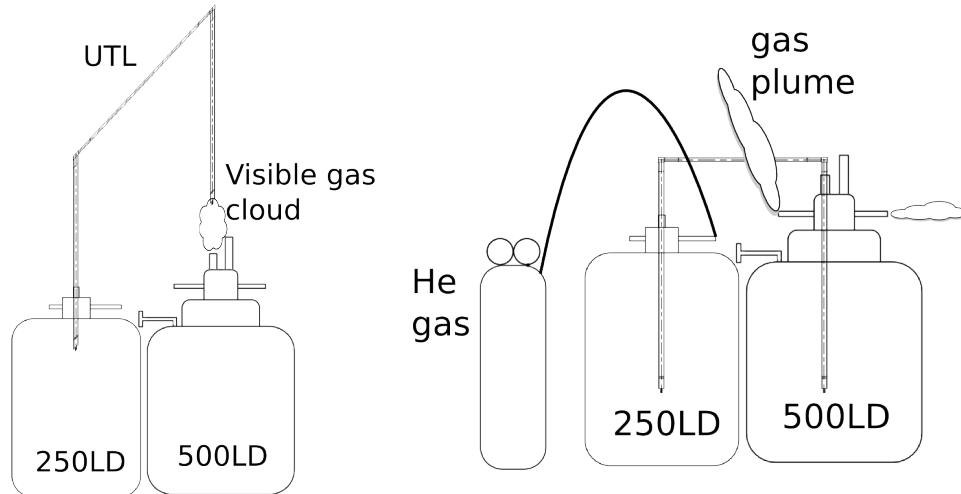


Figure 4.8: Inserting the UTL first in the 250LD, then the 500LD, as described in Steps 7 and 8. The pressurization line is attached to the 250LD to facilitate liquid transfer.

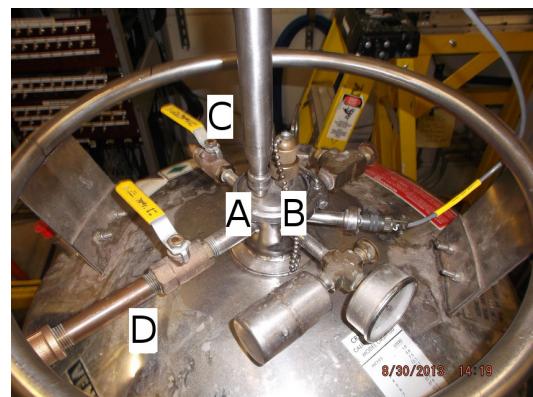


Figure 4.9: Photo of 100LD: A) LHe inlet port, B) LHe outlet port, C) pressure relief valve, D) pressurization/exhaust port.

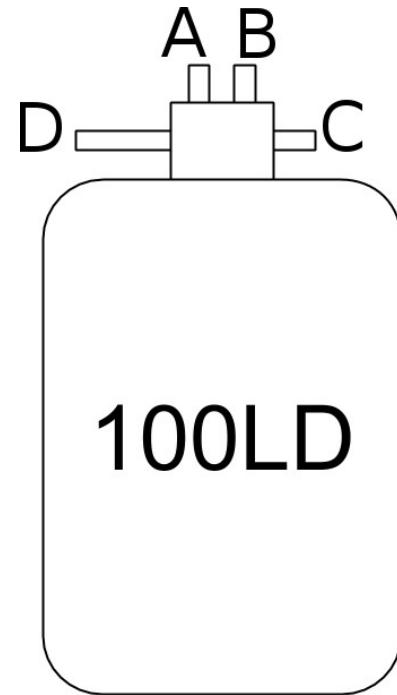


Figure 4.10: 100LD cartoon.

1. Purge the 500TL with helium gas. Make sure the flow valve is open.
2. Make sure the correct Goddard fittings are on the 500LD and/or 500TL.
3. Install the 100TL in the 100LD and attach the bayonet to the WTL.
4. Close the magnet LHe inlet valve on the WTL (blue valve).
5. Open the 100LD exhaust valve, close the pressure release valve and bung the main flow valve.
6. Remove the bung from the WTL by the pump station area.
7. Hook up the helium gas line in the GV to the 100LD exhaust port and pressurize to 2 PSI.
8. Verify there is a positive pressure on the WTL by the pump station. Increase the helium cylinder pressure if necessary. Wait 10 minutes for the 100LD, 100TL and WTL to purge.
9. Close the 500TL flow valve and stop purging it.
10. Vent the 500LD through the pressurization port.
11. Open the 500LD main flow port.
12. Close the 500LD pressurization and pressure relief valves.
13. Slowly lower the 500TL about 30 cm until helium gas starts flowing out. It may need to go further in first, depending on how full the 500LD is.
14. Watch the gas plume coming out of the end of the 500TL. Slowly lower the 500TL in the 500LD at a rate of 1cm/s until the plume looks like a blowtorch, then immediately insert into the WTL and tighten.
15. **Immediately remove the pressurization line from the 100LD exhaust port and open the 100LD pressure relief valve.** Verify gas is coming out the exhaust and make sure the exhaust is pointing towards the center of the room and not at personnel (the plume will quickly grow).
16. Slowly lower the 500TL all the way down into the 500LD and tighten the Goddard fitting. Use the heat gun if the Goddard fitting freezes before the 500TL is completely in.

17. Hook up the pressurization line to the 500LD pressurization port using the procedure in Section 5.14 and open the 500LD pressurization valve.
18. Turn on the 100LD level probe and watch the monitor for the level to increase.
19. When the 100LD is full there will be a large white plume coming from the exhaust and the level probe will read around 24.9 (arbitrary units). Close the 500LD pressurization valve and remove the pressurization line. Then vent the 500LD through the pressurization port and open the 500LD pressure relief valve.
20. Close the 500TL flow valve.
21. Close the 100LD exhaust valve.

100LD Fill During Run

1. Verify the 500LD is full enough to top off the 100LD. If not see above 500LD fill procedure.
2. Set the regulator on the helium pressurization cylinder to 2 PSI.
3. Open the 100LD exhaust valve and make sure the exhaust port is safely facing the center of the room.
4. Hook up the pressurization line to the 500LD pressurization port using the procedure in Section 5.14.
5. Open the main flow valve on the 500TL, then close the 500LD pressure relief valve.
6. Open the 500LD pressurization valve.
7. Verify the 100LD Goddard Fittings are not leaking gas. If they do, remove the pressurization from the 500LD and close the 500TL main flow valve immediately. Consult a high ranking lab official.
8. When the 100LD is full, remove the pressurization from the 500LD, open the 500LD pressure relief valve and close the 500TL main flow valve.

4.2.2 Cool ^4He Section

When the 100LD is full, the fridge has been backfilled with helium gas and the rest of the system is ready for a cooldown, the LTL is inserted into the 100LD and fridge. The separator pumps on the fridge are started and the pressure difference is enough to pull helium from the 100LD (no dewar pressurization needed).

The LHe travels through the LTL into the separator, which collects liquid while rejecting gas. The gas is pulled out through two lines, the “separator” line which precools the ^3He inlet, and the shield line which provides a heat sink for the radiation shield.

When the separator is full, control valve 1 (CV1) allows helium to flow to the evaporator and control valve 2 (CV2) directs LHe to the holding magnet and IVC can. CV2 is also used to quickly cool the MC when target material is at risk of thermal damaged.

The evaporator can hold LHe at 4K, but evaporator pumps are run to lower the helium vapor pressure and lower the temperature down to a minimum of about 1K.

Insert LTL

1. See Section 4.1.1 to make sure the fridge is ready to take LHe.
2. Check the 100LD level probe and top off if necessary (see directions above).
3. Blow out LTL with helium gas.
4. Make sure 100LD can roll out from under the fridge stand area (so the LTL can be inserted), back under the fridge stand area, and then up to the fridge. The IVC pump cart may have to be shifted around.
5. Put 4 PSI on the fridge via the following steps: adjust the helium gas regulator to 4 PSI, open the valves between the regulator and the separator manifold, close all other separator manifold valves (so gas does not flow to separator pumps).
6. Remove the LHe inlet bung and KF clamp from the fridge and verify the positive pressure from the separator.
7. Stop purging LTL with gas and place a rubber in the fridge side.

8. Remove bung from 100LD main flow port and slowly insert the LTL until it lines up with the fridge.
9. Roll the 100LD in place under the fridge stand.
10. Remove the LTL bung and wait for the gas plume to look like a blue torch flame, then roll the 100LD towards the fridge until the LTL is fully inserted and tighten the KF clamp around the LTL-fridge inlet joint.
 1. Open CV1 and CV2 about 2 turns each.
 2. Immediately turn off the pressurization from the helium cylinder and turn on the separator pumps. Open all valves between the separator pumps and the fridge, including bypass valves. Pump on shield at max flow for 30 minutes and separator for 15 minutes[Nii13].
 3. Separator flow is only kept high when there is no liquid in the separator yet; thereafter control the shield flow keeping it below 20 mmol/s and reduce the separator flow to 5 mmol/s (see Appendix B for flow conversion).
 4. TODO Cooldown diverges from June 2014 cooldown at this point; consider separate document

4.2.3 Cool ^3He Section

4.2.4 Dilution

4.3 Polarization

4.4 NMR

4.5 Frozen Spin

Loreum ipsum dolor sit amet, consectetuer adipiscing elit. Etiam lobortis facilisis sem. Nullam nec mi et neque pharetra sollicitudin. Praesent imperdiet mi nec ante. Donec ullamcorper, felis non sodales commodo, lectus velit ultrices augue, a dignissim nibh lectus placerat pede. Vivamus nunc nunc, molestie ut, ultricies vel, semper in, velit. Ut porttitor. Praesent in sapien. Loreum ipsum dolor sit amet, consectetuer adipiscing elit. Duis fringilla tristique neque. Sed interdum libero ut metus. Pellentesque

placerat. Nam rutrum augue a leo. Morbi sed elit sit amet ante lobortis sollicitudin. Praesent blandit blandit mauris. Praesent lectus tellus, aliquet aliquam, luctus a, egestas a, turpis. Mauris lacinia lorem sit amet ipsum. Nunc quis urna dictum turpis accumsan semper.

Chapter 5

Procedures

5.1 Evacuating Transfer lines

Helium transfer lines (TLs) are vacuum insulated tubes. One end is plunged directly into a helium dewar while the other delivers liquid helium to a vessel, with a pressure differential driving the helium through. If the vacuum jacket is compromised, the helium will cool the outer shell of the TL and frost up (helium will of course stop transferring). Once this happens, the only solution is to remove and warm up the transfer line and evacuate the vacuum space.

TLs are evacuated with a leak checker until the leak rate plateaus. It is recommended to pump overnight, but in a bind it can be done quicker. An evacuated TL will be fine for months, but in practice it is safer to pump down each TL a few days before beginning a cooldown in case the vacuum space was somehow compromised in storage/transport.

5.1.1 Cryofab TLs

Cryofab TLs have a proprietary vacuum pump out ports which require a special adapter to KF (Figure 5.1).

To pump out the vacuum space:

1. Warm up the LD.
2. Set the TL on the ground.
3. Loosen the brass fitting on the adapter (do not loosen all the way or the o-ring will fall out).



Figure 5.1: The Cryofab transfer line vacuum space to KF-25 adapter with stem pushed all the way in. The brass nut fits a 1-3/16" wrench.

4. Pull the adapter stem out about an inch and fit the adapter on the TL pump out until almost the entire vacuum jacket pump out is covered (Figure 5.3).
5. Tighten the brass nut.
6. Support the weight of both ends of the TL on stable objects and attach the KF adapter to the leak checker (we use a chair and lab bench in Figure 5.4). Do not allow the TL to torque about the adapter or the TL and adapter may both be damaged.
7. Push the stem in.
8. Screw in the stem about 5 turns (count them if this is your first time) while gently pushing in so the threads grab. The threads are tightening when the stem is pulled in slightly each turn.
9. Gently pull the stem out half an inch. **If you pull too far you will break the adapter; see Figure 5.5.** If TL was under vacuum, you will hear a hiss indicating the vacuum jacket is vented.
10. Start the LD pump and wait for 2 minutes. This should bring the LR to a low enough value to detect an external leak. Spray both vacuum connections on the adapter with helium.
11. If the KF joint leaks, clean it and try again. If the TL-to-adapter joint leaks, either a) the adapter is not entirely over the TL pump out port or b) the brass nut is not tight enough. Try again.
12. Pump until the LR stops declining, preferably overnight.

13. Again, spray both joints with helium gas. Start over and remake the offending connection if necessary.
14. Flow helium gas through the TL (if there is a valve like Figure 5.2, open it first). If there is a leak between the inner cavity of the TL and the vacuum space the TL is broken and not suitable for a cooldown.

At the end of leak checking, use the following procedure to disconnect the TL from the LD:

1. Push the stem straight in all way.
2. Vent the LD. If there is a long hiss, the stem was not pushed in all the way and the vacuum space is now filled with air. A short hiss indicates only the small volume in the KF adapter was vented (this is good).
3. Unscrew the stem. If you counted the number of turns in the above section, you know how many times to unscrew it. Otherwise, just turn a few dozen times to be sure. Not entirely unscrewing the stem will result in the TL being entirely vented in the next step.
4. Pull the stem back about an inch. **If you pull too far you will break the adapter; see Figure 5.5.** You should not hear a hiss.
5. Disconnect the KF connection from the LD and set the TL on the ground. This may require two people to do safely depending on the weight of the TL. Be careful not to bump the TL.
6. Loosen the brass nut about two turns and wiggle the adapter off the TL. Unscrewing the nut too much will make the connector fall apart; **unscrewing too little will break the adapter and/or TL if you try forcing it off.** Do not try forcing it.

5.1.2 LTL

The LTL, named for its shape (Figure 5.7), was designed at CERN. It did not come with a vacuum pump out adapter, so we made one by modeling the brass screw that plugs the pump port.



Figure 5.2: A: to TL stinger; B: to TL bayonet; C: emergency popoff valve; D: LHe flow valve; E: vacuum jacket pump out.



Figure 5.3: Adapter on vacuum jacket pump out.



Figure 5.4: Stem pushed in.



Figure 5.5: Stem pulled out.



Figure 5.6: The LTL adapter.



Figure 5.7: The LTL.

Gather the following:

- LTL
- KF adapter to LTL
- teflon tape
- flathead screwdriver
- leak detector (LD) with KF25 oring/clamp
- ziplock bag

1. Warm up the LD.
2. Wrap the KF adapter's threads with one or two layers of teflon tape.
Do not block the tapered end with tape or it might ruin the seal.
3. Remove the top screw from the LTL and place it in the ziplock bag.
4. Open side screw valve 3 turns to break the vacuum. If there is a hiss,
the LTL was holding a vacuum.
5. Hand tighten the adapter in the top hole while being careful not to
let the emergency popoff valve fall out. It might be tricky to get the
threads right with the tape between them, but do not strip the threads!

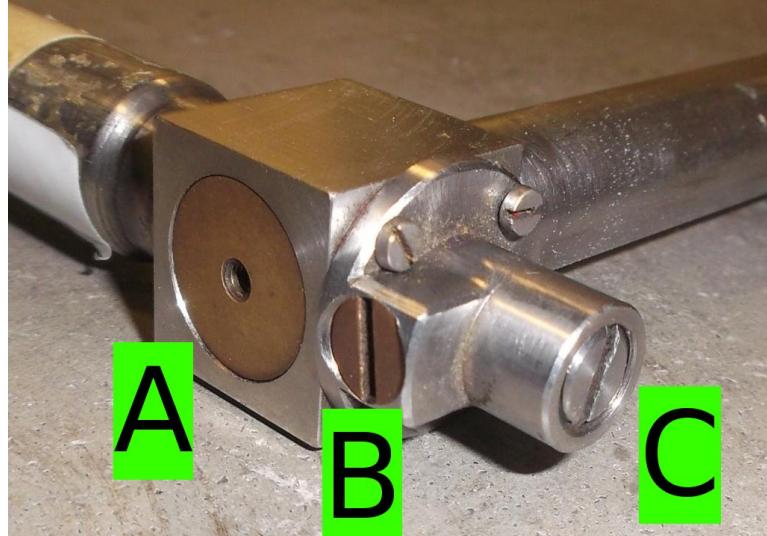


Figure 5.8: The three ports on the LTL: A) emergency relief popoff (this will shoot very high when the LTL freezes!), B) vacuum jacket pump port and brass placeholder screw, C) side screw valve.

6. Start pumping with the LD. The leak rate should approach 10E-8.
7. After a few minutes of pumping, spray the KF and threaded joints with He gas to make sure the connection is leak tight. If the threaded joint is leaking, **carefully** rotate the LTL $1/8$ of a turn, up to $1/4$ turn max. If the joint is still leaking, the teflon tape will have to be removed and reapplied.
8. Pump overnight.
9. Close side screw valve and tighten.
10. Stop and vent LD and break the KF connection.
11. Unscrew the LTL adapter, making sure there is no teflon left over in the LTL port, and screw the brass top screw back in place.
12. Remove all teflon tape from the adapter.



Figure 5.9: Pumping down the LTL. The side screw valve is opened to allow the LD to pump the vacuum jacket.

5.2 Indium Extrusion

Hifrost requires indium wire precisely 1 mm in diameter for making MC and IVC seals. Indium seals must be scraped away after use, and the metal is always recovered (due to its scarcity). The Indium Association of America charges roughly \$1000 for 13 feet, or about \$20-30 per attempted seal, so it is worth creating our own wire from the so-called “scraps”.

Indium extrusion is the name of the process that transforms bulk indium into nicely shaped wire.

5.2.1 Background

Indium scraps or ingots are put into a hollow cylinder with a small hole in the bottom, and the indium is softened with heat and pushed through to make a 1 mm wire. The melting point of indium is about 150° C, but approaching this temperature leaves the indium “runny” and impossible to shape. We have found heating the cylinder to 100° C with a heating tape yields the ideal consistency.

Our first indium die, the plate at the bottom of the hollow cylinder that has a small hole for shaping indium, had a 1 mm hole. For whatever

reason, this produced indium 0.75 mm in diameter (too small for the HiFrost flanges). We increased the die to 41 mils (about 5% larger), and the indium wire came out at 1 mm.

When extruding the indium, a varying force on the hydraulic pump sometimes leads to uneven wire (something like sausage links). To prevent this, a continuous, non-varying force should be applied for an entire press.

5.2.2 Materials

The **extrusion cylinder** is an aluminum cylinder with six 1/4-20 bolt holes tapped in the bottom. The inside of the aluminum is lined with another hollow cylinder made of steel, which precisely matches a **steel rod** that slides in and out. If the rod and inner cylinder do not precisely fit, either the rod will not fit in the cylinder or the indium will leak out during the press. There is a small hole drilled in the side of the aluminum for a thermocouple.

Screwed into the bottom of the extrusion cylinder is the **indium die**, a plate with a hole through the middle to push indium through. The entrance hole in top of the die matches in the extrusion cylinder, and the exit hole in the bottom of the die shapes the indium wire. The die is a separate component from the extrusion cylinder for cleaning and modularity reasons.

Indium is generally recovered from scraps of old seals or ingots are bought from the Indium Association of America (currently a pound costs \$300). Ingots should be cut with a knife into smaller blocks to fit inside the extrusion cylinder.

The **extrusion stand** is a mounted, horizontal bar for the hydraulic press to push against so pressure is applied to the extrusion cylinder. There must also be a hole in the base of the stand for the pressed indium wire to be collected.

A **thermocouple** sensor is used with an accompanying **thermocouple readout** to measure the temperature of the extrusion cylinder.

Heat tapes warm the indium so it may be pressed. The tapes should be long enough to wrap around the cylinder to warm it evenly.

A layer of **aluminum foil** is wrapped around the heating tapes to keep the heat from dissipating away.

Latex and cloth gloves are for personnel safety and keeping the indium sanitary.

A **variable alternating current** unit adjusts the temperature of the heating tape.

The indium is extruded onto sterile **Absorbond Wiper Sheets, TXTM** 9"x 9", 409, and kept for storage.

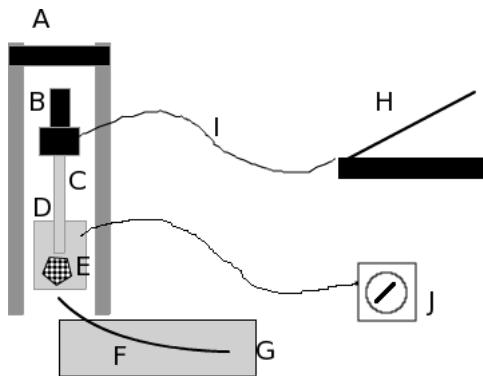


Figure 5.10: Indium extruder setup: A) extruder stand, B) hydraulic pump, C) steel rod, D) extrusion cylinder, E) indium scrap/ingot fragments, F) finished wire exiting die, G) sterile wipe sheet, H) hand jack, I) hydraulic line, J) thermocouple gauge

An *EnerpcTM* P39 Porta Power **hydraulic hand pump** provides the force for extrusion.

5.2.3 Extrusion process

1. Prepare connections for heating tape and thermocouple gauge and readout.
2. Remove steel rod from extrusion cylinder and insert blocks of indium.¹
Place the steel rod back in the cylinder.
3. Wrap the cylinder with heating tape and aluminum foil.
4. Place the cylinder/rod under the extrusion stand.
5. Set the variable current source so the heating tape reaches about 70° C, which takes around 20 minutes. Increase the temperature in steps, allowing the thermocouple to equilibrate between steps, to 100° C.
NB: on our VariAC unit the 50 mark corresponds to 75° C and the 60 mark corresponds to 100° C
6. Set the hydraulic jack on the extrusion cylinder and under the stand's horizontal bar.

¹After losses, 100 grams of indium provides approximately 10 feet of wire instead of the theoretical 50 feet.

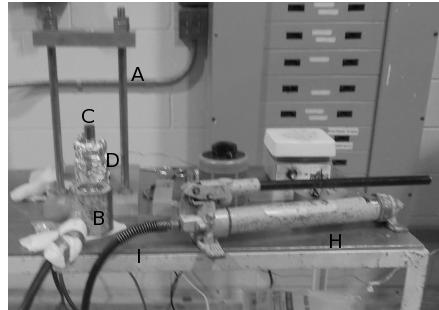


Figure 5.11: A photo of the setup.



Figure 5.12: The extruder cylinder without the die. The steel rod is plunged all the way down and leftover indium can be seen on the surface.

7. Use the hand jack to slowly push the indium through the extrusion cylinder and out the die. The press takes about 5 minutes. A second person slowly winds the indium wire on the wiper sheets.

5.3 Making Indium Seals

1. First, clean both flanges the indium will seal together.
2. Measure the diameter of the indium to be used; it should be 1 mm
3. Cut a length just greater than the circumference of the flange (MC: 158 mm; IVC: 271 mm) plus a centimeter or two. Do not allow the indium wire to touch any dirty surfaces.
4. Wrap the wire around the proper aluminum setting tool.
5. Using a razor blade, cut the wire at an angle as in Figure 5.14, firmly placed on the aluminum block.
6. Flip the mating can (IVC or MC bell housing) upside down and press firmly to the aluminum block TODO photo.
7. Turn the parts so the aluminum block is on top, and carefully remove the block so the indium is transferred to the can. Check to make sure the wire was not shifted, especially where the two ends meets.



Figure 5.13: The hydraulic press sitting on top of the extrusion cylinder.

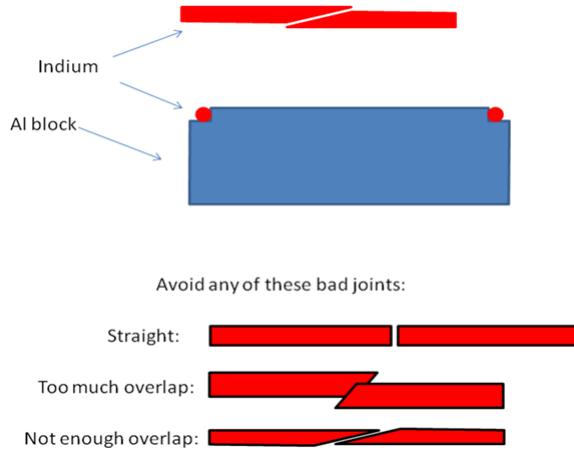


Figure 5.14: How to cut indium wire for fridge seals. [Zha13]

8. When it is time to make the seal on the fridge, carefully raise the can up to the fridge within 1 cm from touching, not allowing anything to touch the indium wire. If the wire is displaced, begin from the beginning. Hold in place with a chemistry clamp. TOOD show photo
9. Insert two appropriate machine screws up through the can, and turn twice so they can support the weight of the can without the chemistry clamp. The indium should not touch the fridge flange.
10. Remove the clamp and turn the remaining four screws 2 times each.
11. Alternate tightening the screws 1/4 turn each in a “star” pattern. That is, do not turn adjacent screws consecutively, and no single screw is tightened twice unless the other five screws have been tightened in between.
12. When the indium is compressed, make sure all the screws are evenly tightened, but do not strip the machine screws’ hex heads.
13. If possible, leak check the indium seal.

5.4 Breaking Indium Seals

1. Remove screws with 1.5 mm allen key and the modified MC wrench. Visually inspect each screw for indium or signs of stripping, and turn

each nut back on the screw it came off of. If there is any torque resistance or sign of stripping, the screw was likely damaged when the indium seal was tightened and both the screw and nut should be discarded. If there is no damage, return the nuts/screws to the bag labeled “MC Screws and Nuts”.

2. Use two MC screws and the 1.5 mm allen key as jacking screws the break the indium seal. Defining the top of the fridge (in the horizontal orientation) as 12:00, the jacking screw holes are at 1:00 and 5:00 look upstream. The jacking screw holes can be identified by threaded holes in the MC bell housing that have no matching through holes on the fridge MC flange. Alternate between screws making 1/4 turns until the MC is free.
3. Remove the MC from the fridge, careful not to damage the sensors on the copper dam.
4. Remove the jacking screws from the MC and place in the bag labeled “MC Screws and Nuts”.
5. Scrape indium, if any, off the MC and place the MC in the safety zone.
6. Scrape indium off the fridge flange, recovering as much as possible in the indium scrap container; be very careful not to touch the stainless steel beam window.

5.5 Weighing LHe Dewars

5.6 Magnet Lead Installation

5.7 Mounting Fridge

5.8 Dismounting Fridge

5.9 Cleaning Viton/Buna O-rings and Grooves

HiFrost only uses isopropyl alcohol for cleaning vacuum parts. Methanol bottles are scattered around DFELL, but they are not to be used.

Most o-ring encountered are for KF joints. It is not practical to clean every o-ring every time it is touched, and this is unnecessary because our

vacuum system is not UHV, so most KF joints can be visually inspected for dirt before use.

All ^3He system o-rings must be leak checked every time they are taken apart, regardless of whether they were cleaned or not.

To clean a KF o-ring:

1. remove o-ring from center piece
2. moisten a Kim Wipe with isopropyl alcohol
3. gently pinch the o-ring with the wipe, rotating it around to clean the whole surface
4. replace the oring over the centerpiece
5. wipe down the outside of the o-ring one more time
6. reassemble the KF joint, trying to touch only the centerpiece as much as possible (tip: gravity is your friend, here)

5.9.1 Fridge flanges

The triple flange, ^3He back flange, evaporator pump out and ^3He pump out all have rubber o-rings set in grooves. If any of these leak the fridge will not be able to achieve dilution (or ^3He will leak out), but they cannot be leak tested during a cold-target run because there is no time to pump down the system with a leak detector.

To minimize the chance of them leaking, the following should be done each time one of these flanges are opened:

1. clean o-ring with isopropyl alcohol as described above
2. dampen qtip with alcohol and swap out o-ring groove (this will be very dirty from vacuum grease)
3. use a minimal amount of vacuum grease to coat the o-ring all the way around
4. replace the o-ring and seal the flange immediately

The vacuum grease feels unpleasant on hands, so latex gloves may be used. Otherwise, be sure to wash up afterward so vacuum grease doesn't contaminate other fridge components.

5.10 Changing Kenol Connectors

5.11 Installing Heating Tapes

5.12 Filling LN₂ Trap

5.13 Check 250LD/500LD LHe Level

The LHe level of helium dewars is checked upon delivery and after they are used to fill the 500LD. The 500LD LHe level is checked after it is filled, after initially filling the 100LD and periodically throughout the run. The 100LD has an internal level probe installed and does not use this procedure.

1. Remove zip ties and plastic wrap on full 250LD.
2. Make sure there is nothing sensitive directly above the dewar, then release pressure by opening main flow valve on top of the dewar. A large, cold plume will exhaust pressurized helium gas.
3. Use cryogenic gloves to slowly lower the level probe into the main flow valve. A smaller cold plume will exhaust upwards as the level probe cools down.
4. Plug the level probe into the level probe readout and record the reading. As a sanity check, lift the level probe a few inches to make sure the readout changes correspondingly.
5. Disconnect the level probe and remove it from the 250LD, storing it where it can safely warm up.
6. Close the main flow port.

5.14 Attach Pressurization Line to LHe Dewars

The 250LD and 500LD have LHe pushed out of them by pressurizing the top of their liquid contents with helium gas. The helium gas pushes down on the liquid surface, pushing liquid up a transfer line and out of the dewar.

The procedure below ensures air will not get in to the pressurization line and freeze the dewar port.

1. Check that the pressurization cylinder is not empty and a full helium cylinder is nearby. Make sure the regulator has single PSI increments discernable.
2. Open the pressurization cylinder and set the regulator to 2 PSI. Verify helium gas is flowing out.
3. Slowly open the dewar pressurization valve on the to-be-pressurized dewar until pressurized gas begins to flow out at about 5 PSI. If the pressurization frosts, close it and wait until it warms up, using a heat gun only if necessary.
4. With gas coming out of both the dewar pressurization port and the pressurization cylinder, mate the two connectors and close the dewar pressurization valve.

5.15 Purging Helium Lines

5.16 Vacuum Rise Testing

A leak checker is the proper way to determine if a vacuum volume is leak tight. However, it is not always practical and sometimes very difficult to leak check a large volume, so we use a vacuum rise test as a coarse characterization of the vacuum fidelity.

The concept is to pull a vacuum down in a volume, isolate it from the pump, then watch if the vacuum rises. Outgassing and a leak are generally the only two effects we will see, but outgassing usually stops the pressure rise far below 1 atm.

To set up, gather the following:

- pump
- associated KF vacuum fittings/hoses
- KF pressure gauge

The procedure is

1. Set up as shown in Figure 5.15 or Figure 5.16. If the system is powered (such as a pump), make sure it is unplugged.
2. Open v and pump down s until the pressure gauge g stops decreasing.

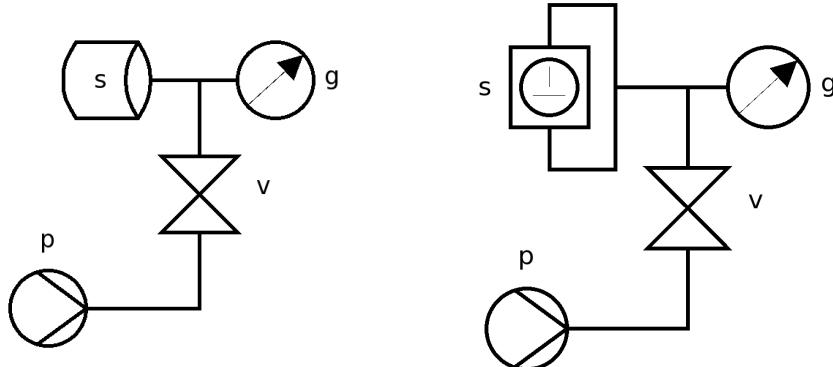


Figure 5.15: Vacuum rise schematic: s is the system to check; g is a pressure gauge; v is a leak-tight manual valve; p exhausts to atmosphere.

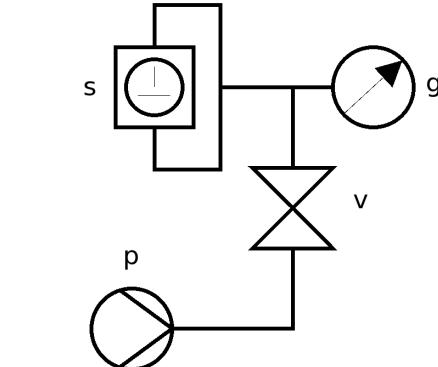


Figure 5.16: The same as Figure 5.15 except s has two ports, the inlet and exhaust of a pump, in this case.

3. Close v and turn off p .
4. Monitor g until it stops rising. This can take from minutes to days.

If the system s is tight it will outgas to tens of mbar (depending on if it is contaminated with oil, water, etc) then stop. If it is necessary (and possible), s should then be checked with a leak detector to verify it is helium tight.

If the pressure gauge reaches an atmosphere, there is a leak in s . The time this takes determines how big the leak is. A quick way to verify s contains the leak is replacing it with a blank off and repeating the vacuum rise test. If g reaches an atmosphere again, there is a leak in the plumbing which should be repaired before trying the vacuum test again.

5.17 Stripping Fridge

5.18 Assembling Fridge

When putting the fridge together, there are two primary concerns: forgetting to install something (e.g., superinsulation, microwave guide support) and leaky seals. This section does not include a cold target load and assumes the fridge is hanging vertically on the fridge stand in the Vault.

5.18.1 MC

Tools

- qtips/wire cutters/isopropyl alcohol
- indium scrap box
- 1.5 mm allen key
- MC open faced allen wrench (4.0)
- MC
- 1 mm indium wire
- chemistry clamp
- 6x M2x12 machine screws and nuts, 316SS

Read the section on making indium seals.

Hold the MC in place with the chemistry clamp so the bell housing is about 1 cm from touching the fridge flange. Place the bolts through the MC holes and gently tighten the nuts one or two turns so the bolts don't fall. Inspect the indium wire to make sure it is still in place, then tighten the nuts in a "star" pattern (never tightening adjacent screws consecutively). Hold the bolts with an allen key while turning the wrench so the heads do not strip. Use only the special MC wrench so the fridge does not scratch.

If loading conditions permit, leak check the MC before moving on.

5.18.2 Microwave Guide Support

Tools

- small flat head screwdriver
- squeeze-to-open tweezers

Three brass screws are usually stored in the threaded struts on the fridge where they fasten the waveguide support in place.

After unscrewing them, raise the microwave guide support over the MC making sure the waveguide is lined up with both the slot in the MC bell-housing and the circular waveguide connector on the fridge. One person holds the support while another tightens the screws. Hold the screws with tweezers to prevent dropping them down towards the MC (if this happens,



Figure 5.17: The MC with a prepared indium seal behind held by a chemistry clamp.

you must start over, careful not to damage the fridge with the loose screw while taking off the support). Warning: the microwave guide support can still fall off when it is being held by one screw.

Make sure the waveguide (or the support) is not touching the MC. Normally this will not happen unless the support has been deformed.

5.18.3 IVC

Tools

- solder station
- 2.5 mm allen key
- indium scrap box
- qtips/wire cutters/isopropyl alcohol
- 6x M3x12 machine screws

Read the section on making indium seals.

Similar to installing the MC, first prepare the IVC indium seal and place it in the groove on the IVC flange. One person holds the IVC within about 1 cm from the fridge while another puts in the screws. Initially, only turn in 2 diametrically opposite screws 1 or 2 turns each so they support the IVC and it no longer needs to be held by another person. Turn in the other 4 screws 1 or 2 turns each. Inspect the indium seal to make sure it is still in place (it should not be touching the fridge flange yet) before turning in the screws in a “star pattern” similar to the MC. Be careful not to strip the screw heads.

Leak check if possible.

5.18.4 OVC

Tools

- pliers for pilot pins
- 5 mm allen key
- three OVC cans (inner, outer, nose cone)
- triple flange orings

- pilot pins
- 2x nut plates
- 8x M6x25 triple flange machine screws

Make sure the orings and grooves are clean (see cleaning orings section). Place the outer OVC can on the blue fridge lift in the Vault, and lower the inner OVC can inside of it, making sure the black lines on the flanges line up. Place the pilot pins up through the pilot holes and secure the cans in place with nuts and threaded rod that fit through the fridge lift platform. Hook up the OVC manifold to the pump out (making sure the hose is coming off at the appropriate angle for Blowfish clearance) and warm up the OVC turbo pump. Raise the nose cone in place, being sure to line up the large dents, and begin pumping. Leak check the OVC.

With the OVC still being pumped, remove it from the fridge lift and raise it over the fridge. Line up the black lines on the triple flange and push the pilot pins in by hand, using pliers if they need to be rotated or pulled out slightly.

Loosely turn in all the triple flange screws, careful not to strip the aluminum taps. Tighten in a star formation so the oring does not become unevenly compressed.

5.19 Removing/Replacing ${}^3\text{He}$ Baffles

Chapter 6

Subsystems

6.1 Thermistor Switch and Display

The thermistor system is set up to provide signals from all of the thermistors to the USB-1608G data acquisition (DAQ) unit (and so be visible and recorded outside the vault), and to maintain the ability to put any of the thermistor readings on a display while in the vault.

The EIO Body thermistor is wired to a dedicated single channel controller, which displays the temperature, has a normally-closed relay with a high-temperature trip point (set to 40 °C right now) in the Cober power supply interlock loop, and provides a DC voltage proportional to temperature to the DAQ (0 to 10 V corresponds to 0 to 100 °C).

Three other thermistors (EIO water flow out, Q-meter P and Q-meter D) are wired to the thermistor switch panels. These are wired in series, and a 100 μ A power supply runs current through these, and DAQ channels are wired in parallel with each of the thermistor wires. The signal is proportional to thermistor resistance. These are 10 k Ω (at 25 °C) NTC thermistors, so the signal is 1 volt at 25 °C. A look-up table will need to be used to convert the readings to temperature.

The switch has 6 positions; Spare 1, EIO Flow, Q-meter P, Q-meter D, Spare 5 and Display Nothing. Set the switch to the thermistor which should be put on the display. When the switch is set to one of the first 5 positions, relays divert that thermistor signal to the display, and connect a short into the 100 μ A PS loop. So, when you look at, say, EIO Flow, the DAQ signal for that channel will be zero volts, and that thermistor's temperature will show on the display. If the switch is set to Display Nothing, all 5 signals go to the DAQ. There is a 29.4 k Ω resistor for the display in the Display

Nothing position, which displays as 0 °C (or maybe 0.1).

The Spare 1 and Spare 5 positions also have 29.4 kΩ resistors, so the display will show zero on those positions as well, and the DAQ will get a pretty constant 2.94 V signal from them. These resistors can be replaced by thermistors in the future should the need for more temperature readings arise.

Referring to the schematic drawing, the switch is 4-Pole 6-Position (4P-6T). Two poles are used to energize the selected relay by sending 24 VDC and 24 Common to the selected coil.

In the not-energized state, each relay connects the thermistor for that channel to the 100 μA power supply. Each thermistor is also wired to a DAQ channel. The thermistors are wired in series, so the same 100 μA flows through all thermistors.

In the energized state, each relay connects the thermistor for that channel to the other two poles of the 4P-6T switch, which connects this thermistor to the thermistor display on the panel. Also, the thermistor is replaced by a short circuit so the 100 μA will continue to flow through the rest of the thermistors.

Also seen is a schematic, Figure 6.1, for the 100 μA power supply, which is based on a μA723 voltage regulator. Operating with a 24 VDC V_{cc} supply, the power supply output can provide 100 μA through a load of 0 to 100 kΩ.

The 6 photos, Figure 6.2, show the front and back of the three small panels which comprise the thermistor system. Switch panel, Display panel, and Relay panel. The Switch panel and Relay panel are connected with a 25 wire DB-25 cable. The display panel also holds a flow switch display which is wired to the EIO cooling water flow switch.

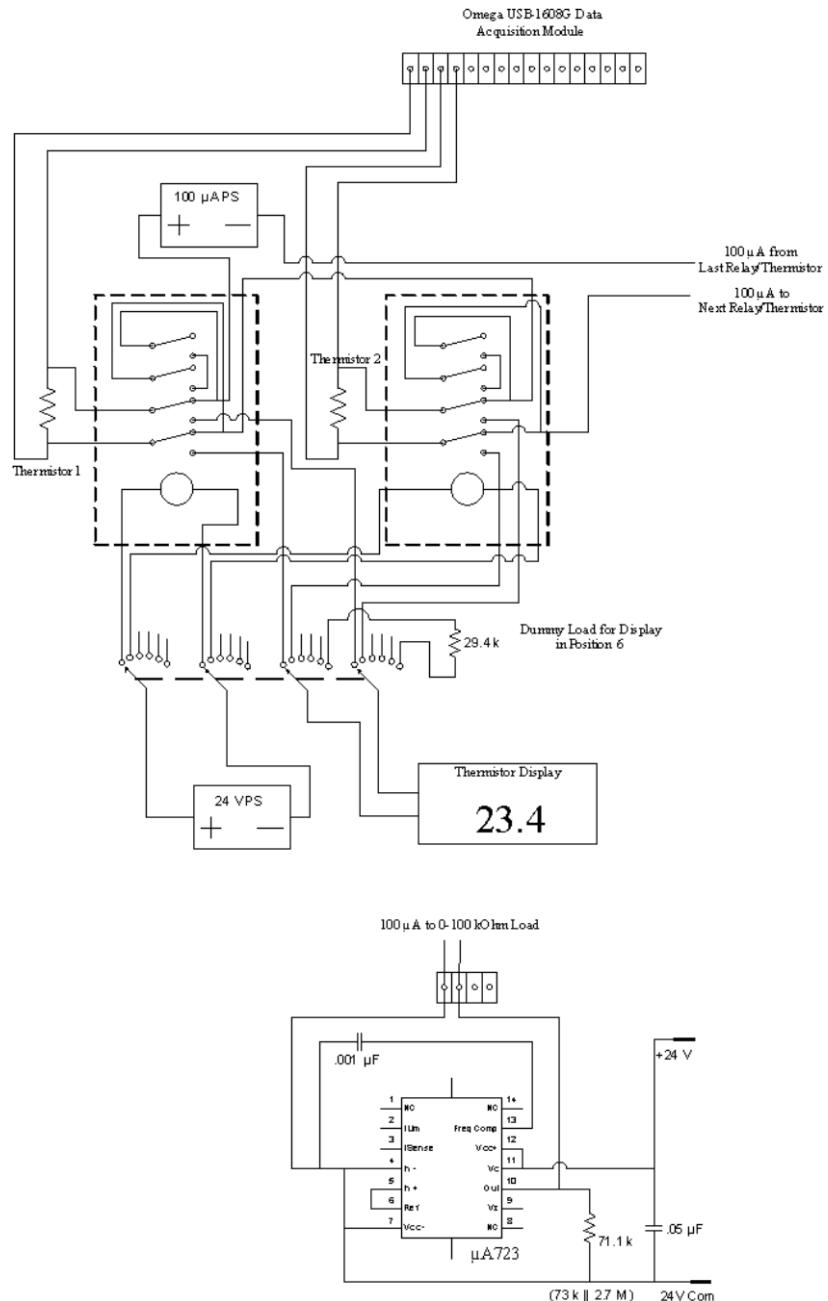


Figure 6.1: Schematic of thermistor system.



Figure 6.2: Photos of thermistor system.

Appendix A

Checklist for cooldown

The preparation for a cooldown is not optional. This checklist should be gone through and double checked before LHe is ordered.

The entire group should be made aware if something on the checklist is not satisfied when LHe is ordered. Subsequent procedures and operations in the Practical Operation (Chapter 4) then need to be reviewed for potential modifications to accommodate the missing prep.

Checklist Before Ordering LHe

Date: _____ Name: _____

Pump station area - ^4He

- flow meters on and working
- laptop webcam working
- $\{\square \text{ UTL}, \square 500\text{TL}, \square \text{ WTL}, \square \text{ MTL}, \square 100\text{TL}, \square \text{ LTL}\}$ pumped down
- $\{\square 250\text{LD}, \square 500\text{LD}\}$ Goddard fittings assembled and ready
- 500LD level probe/monitor plugged in and powered on
- helium gas cylinder hooked up to 500LD
- at least one full backup helium cylinder nearby
- dewar scale and weight record form is ready
- system is purged and backfilled with helium gas

Pump station area - ^3He

- LN₂ trap regenerated
- LN₂ trap filled
- vacuum test ^3He lines
- record ^3He tank pressures in log book (T_1, T_2, T_3)
- flow meter on ^3He gas rack powered on and recording
- system is purged and backfilled with helium gas

Gamma Vault

- 100LD Goddard fittings assembled and installed on 100LD
- $\{\square 100\text{LD LHe}, \square \text{ magnet LN}_2, \square \text{ magnet LHe}\}$ level sensors reading
- fill LN₂ in magnet space
- OVC assembled and pumped down to $\sim 1\text{E}-5$ mbar
- IVC manifold assembled and hooked up to fridge
- Lakeshore sensors read to computer
- AVS-47 bridge is warmed up and all sensors working
- walkie talkies charged

Gamma Vault - Microwaves

- EIO water is running
- microwaves tested with power meter
- wave guide is hooked up to fridge
- thermistors and water flow readouts work and read to computer

Gamma Vault - NMR

- signal achieved with oscillator crystal in mixing chamber coil
- NMR water is running
- PDP running

Appendix B

SLPM Conversion

Some HiFrost documents from CERN refer to flow rates of ^4He and ^3He in millimoles per second. The formula to convert this flow rate to SLPM is

$$[\text{SLPM}] = \left(\frac{60 \text{ s}}{\text{min}} \right) \left(\frac{x \text{ mol}}{\text{s}} \right) \left(\frac{4 \text{ g}}{\text{mol}} \right) \left(\frac{1 \text{ L}}{0.1785 \text{ g}} \right) \quad (\text{B.1})$$

for ^4He and

$$[\text{SLPM}] = \left(\frac{60 \text{ s}}{\text{min}} \right) \left(\frac{x \text{ mol}}{\text{s}} \right) \left(\frac{3.01 \text{ g}}{\text{mol}} \right) \left(\frac{1 \text{ L}}{0.135 \text{ g}} \right) \quad (\text{B.2})$$

for $^3\text{He}[\text{Lin}]$, where $[\text{SLPM}]$ is the standard flow volume, x is the flow rate in mol/s, k is the Boltzmann constant and the definitions for STP are a temperature of 273.15 K and a pressure of 100 kPa.

Figure B.1 shows some values of SLPM to mmol/s flow.

SLPM	mmol/s (³ He)	mmol/s (⁴ He)
5	3.74	3.71
10	7.48	7.44
15	11.21	11.16
20	14.95	14.86
25	18.69	18.60
30	22.42	22.31
40	29.90	29.75
50	37.38	37.19
60	44.85	44.63
80	59.80	59.50
100	74.75	74.38

Figure B.1: Flow rate conversion between mmol/s to SLPM for helium.

Appendix C

Transfer Line Drawings

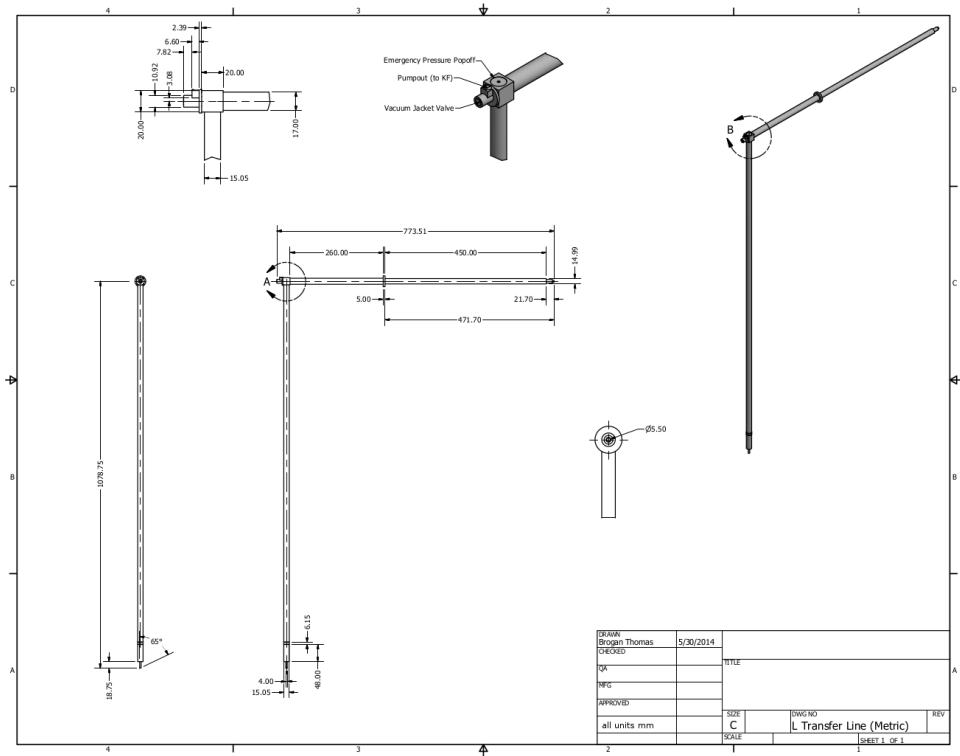


Figure C.1: Drawing of the LTL.

APPENDIX C. TRANSFER LINE DRAWINGS

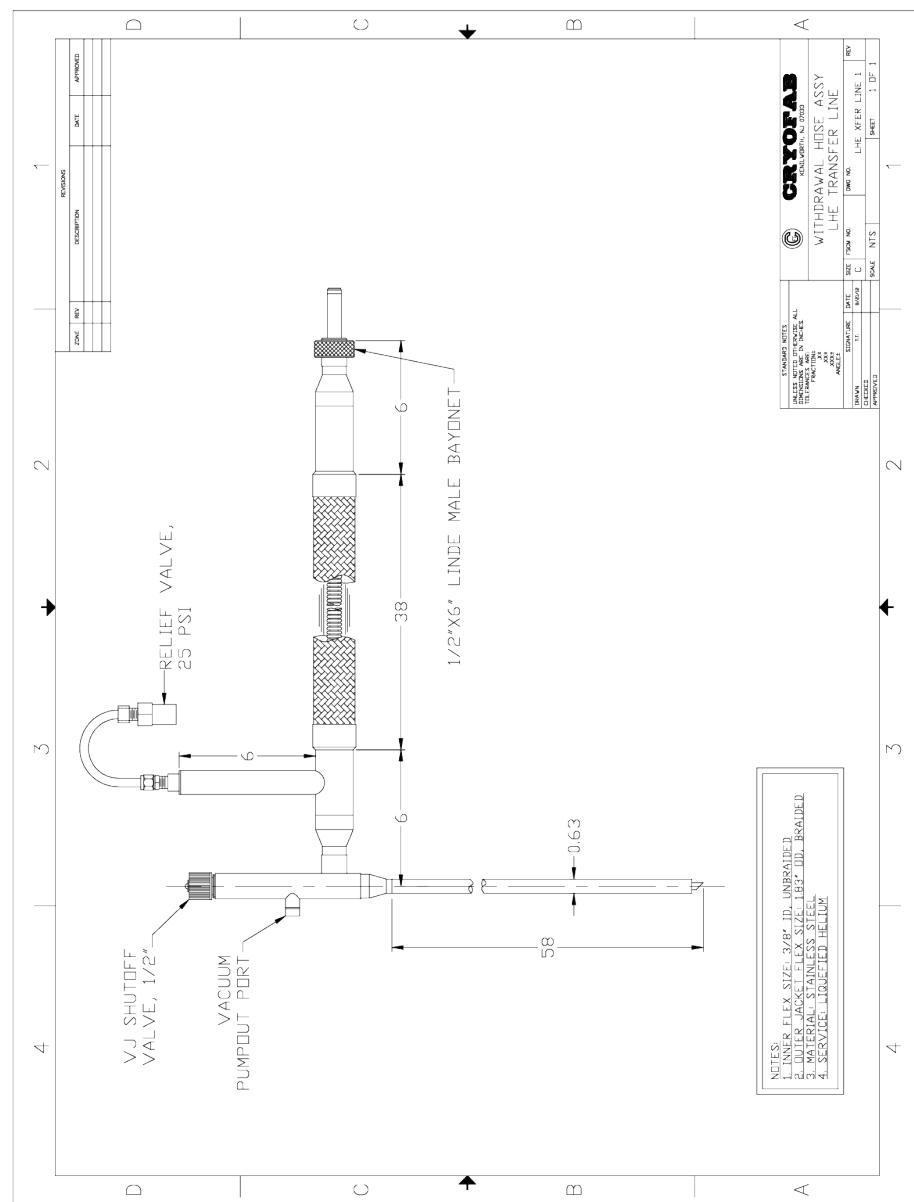


Figure C.2: Drawing of the 500TL.

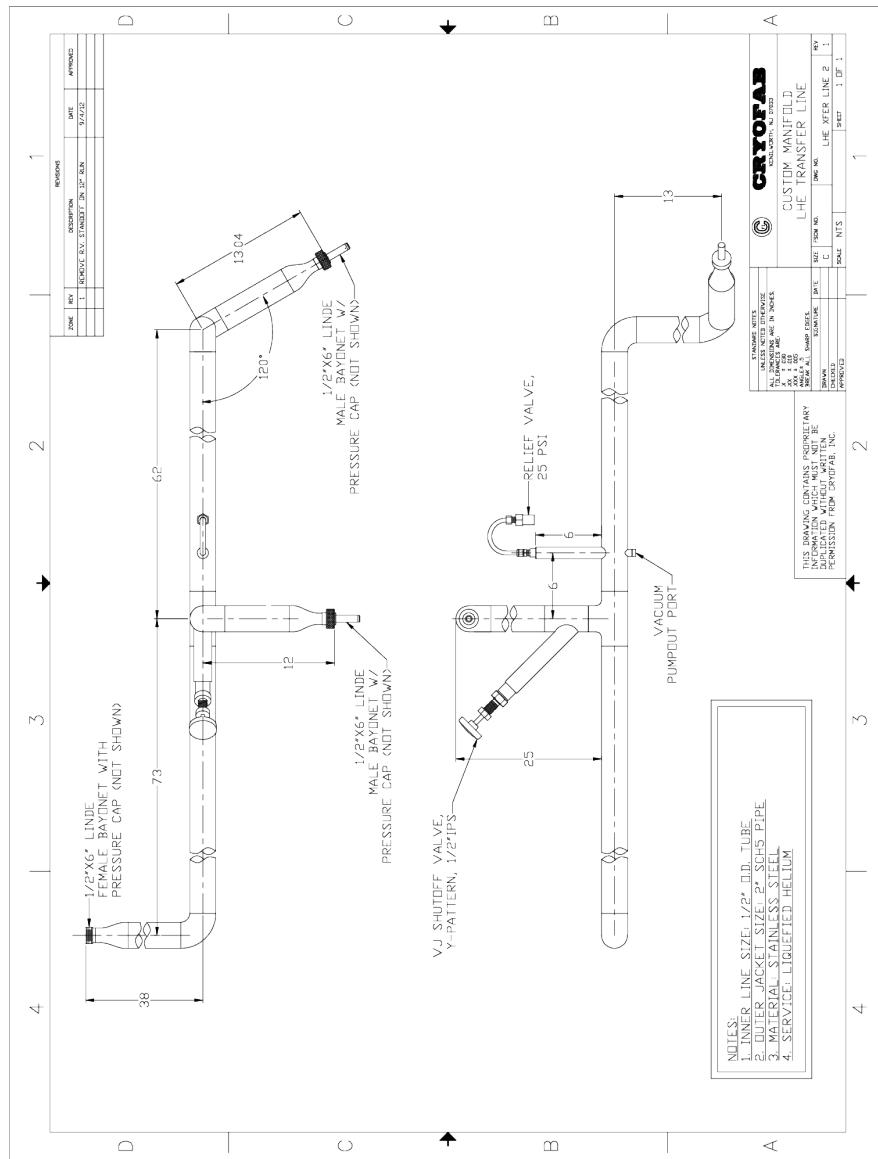


Figure C.3: Drawing of the WTL.

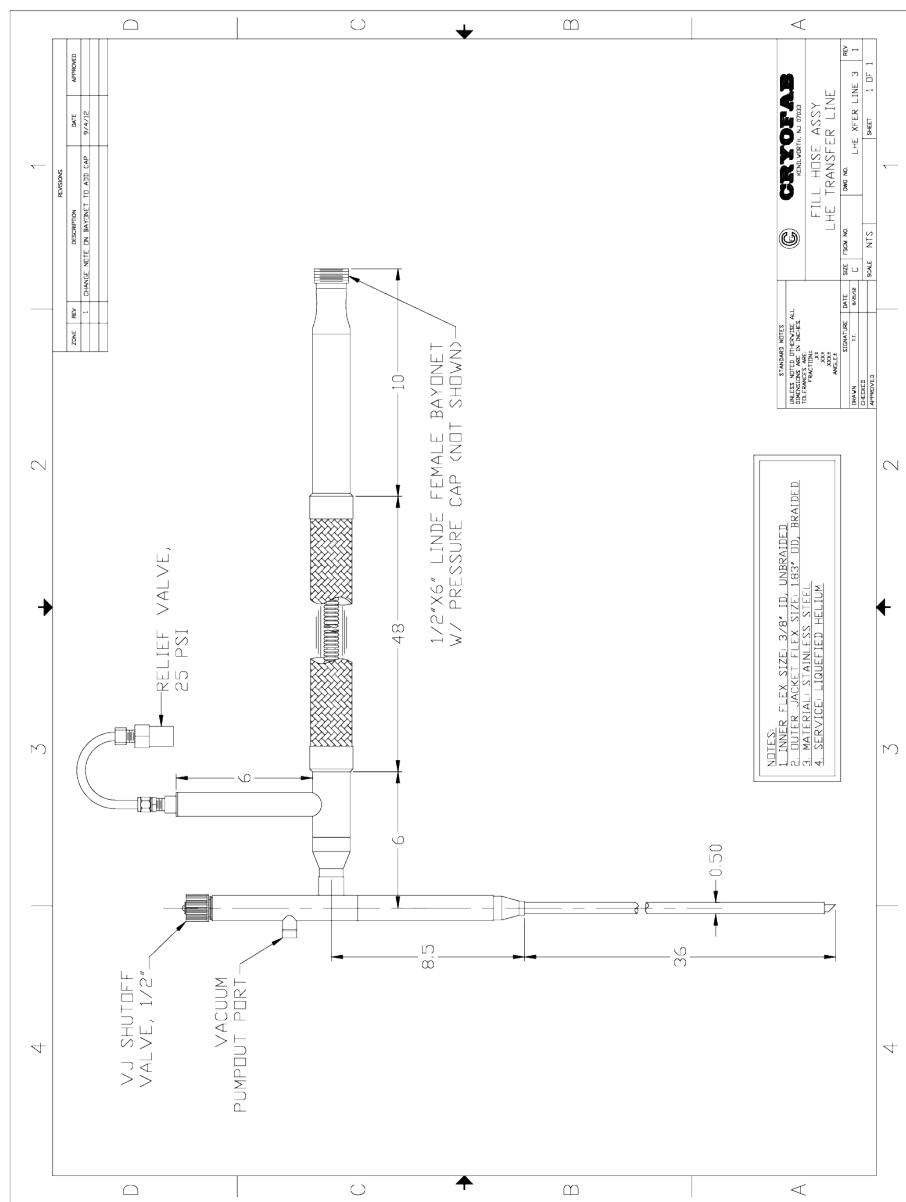


Figure C.4: Drawing of the MTL.

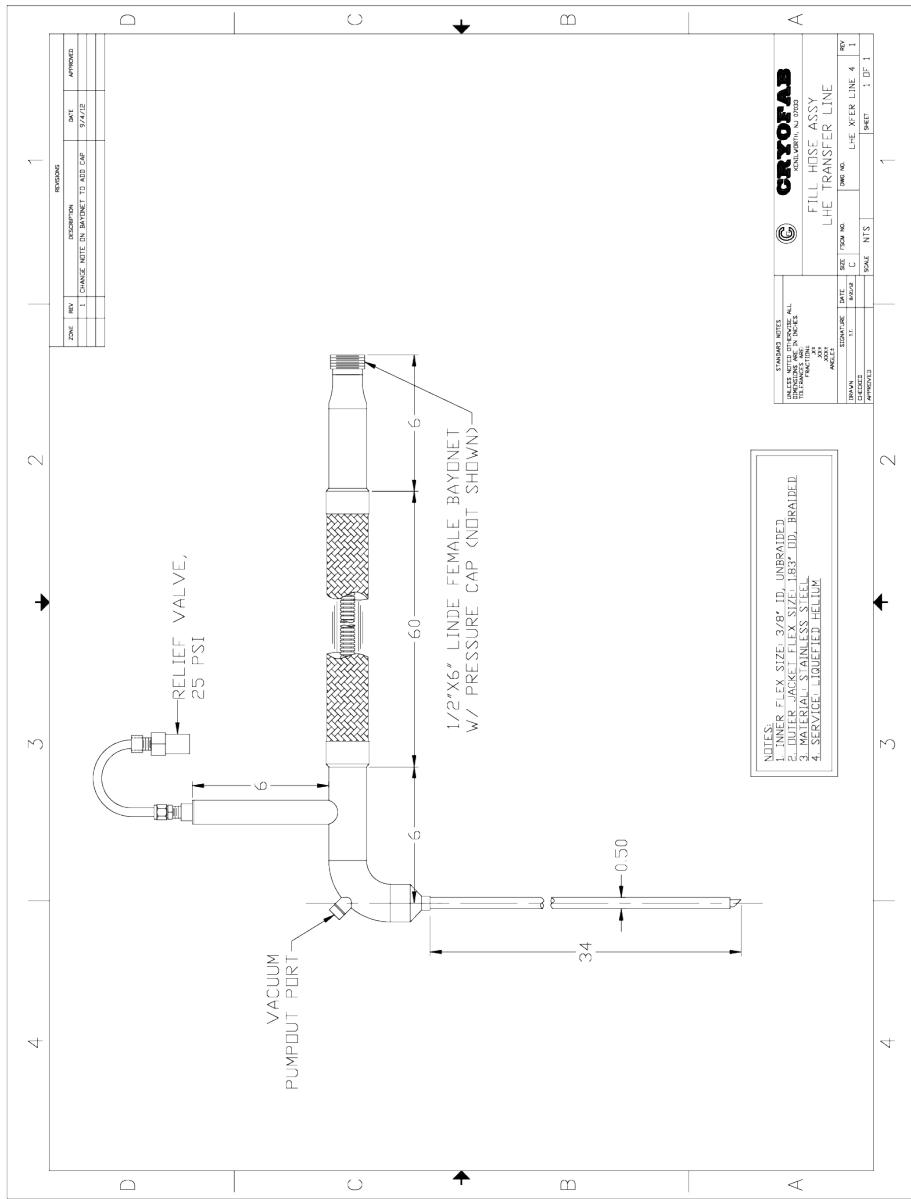


Figure C.5: Drawing of the 100TL.

Bibliography

- [Air] Airgas, *Isopropanol material safety data sheet*, <http://www.airgas.com/documents/pdf/001105.pdf>.
- [Bet89] DS Betts, *An introduction to millikelvin technology*, Cambridge University Press, 1989.
- [CG80] G Court and DW Gifford, *Liverpool nmr module, q-meter fundamentals, and circuit testing*, Ultra Physics LTD, November 1980.
- [CM97] D Crabb and W Meyer, *Solid polarized targets for nuclear and particle physics experiments*, Annual Reviews Nuclear and Particle Physics **47** (1997), 67–109.
- [EPA13] United States Environmental Protection Agency EPA, *Learn about lead*, <http://www2.epa.gov/lead/learn-about-lead>, July 2013.
- [Epl14] EOHW Director Carol Epling, *Cryogenic burn response*, private email correspondence, April 2014.
- [Hoc] Kai Hock, *Cooling techniques*, hep.ph.liv.ac.uk/~hock/Teaching/2011-2012/4-cooling-techniques.pdf.
- [IC] Duke-UNC Brain Imaging and Analysis Center, *MRI Safety Tutorial*, www.biac.duke.edu/research/safety/tutorial.asp.
- [Lin] Linde Electronics and Specialty Gases, *Helium-3 chemical specification sheet*, lindepremiumproducts.com/internet.lg.lsg.usa/en/images/he3-helium-3892_88295.pdf.
- [Mat06] Brent Mattox, *Investigative report on chemistry 301a cylinder explosion*, Tech. report, Texas AM University, January 2006.
- [Nii13] Tapio Niinikoski, *Cooldown with LHe*, private email correspondence, November 2013.

- [OO] Duke Occupational and Environmental Safety Office, *Working Safely with Stepladders at Duke*, www.safety.duke.edu/ohs/Documents/Laddersafety.pdf.
- [RWDL] S. Ryf, T. Wolber, F. Duru, and R. Luechinger, *Interference of neodymium magnets with cardiac pacemakers and implantable cardioverter-defibrillators*, <http://www.ncbi.nlm.nih.gov/pubmed/18334784>.
- [SA] Occupational Safety and Health Administration, *Stairways and Ladders*, www.osha.gov/Publications/osha3124.pdf.
- [Sta98] Kevin C. Staudinger, *Occupational lead poisoning*, <http://www.aafp.org/afp/1998/0215/p719.html>, February 1998.
- [TJNL] Industrial Hygiene Group Thomas Jefferson National Laboratory, *Oxygen deficiency hazard safety booklet*, jlab.org/accel/safety1b/ODH-book.pdf.
- [Zha13] Jinkui Zhao, *On making indium seals*, private email correspondence, June 2013.